




Quartet™ TA Quick Reference Guide

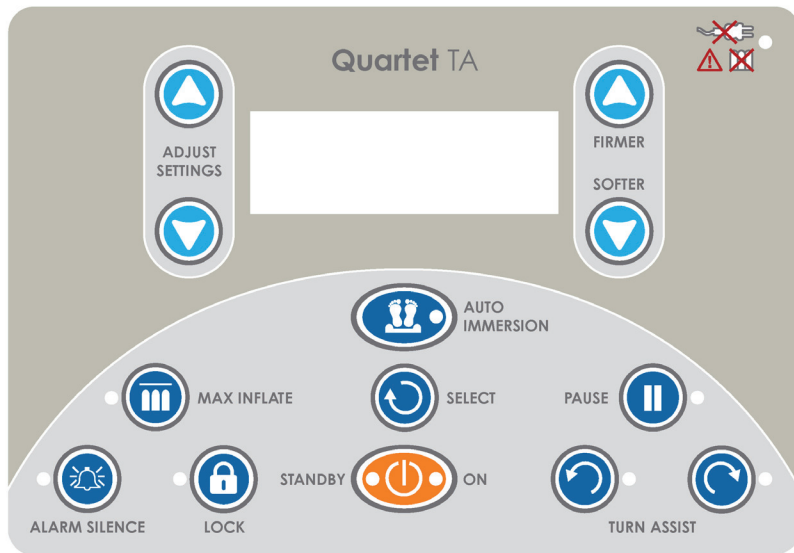

POWER

Press the POWER button. The STANDBY (amber) LED will turn off and the ON (green) LED will turn on.



TURN ASSIST

Press and hold one of the TURN ASSIST buttons and use the SELECT button to toggle between options:



- Centering Function
- Turn Degree
- Turn Time

Use the ADJUST SETTINGS buttons to set the above options. To exit the menu, press the PAUSE button. Press the same TURN ASSIST buttons and the surface will begin turning.



PAUSE

Press the PAUSE button to halt mattress rotation. To set the pause time, hold the pause button for 2-3 seconds and use the ADJUST SETTINGS arrows to set the the pause time between 5-120 minutes (default is 15 minutes). To exit the menu, press the SELECT button.



AUTO IMMERSION

With patient on the mattress, press the AUTO IMMERSION button and the blue LED will turn on. The surface will monitor and automatically adjust immersion settings to match patient weight and position.



LOCK

Press and hold the LOCK button until the LED light turns on, to prevent tampering of settings. To unlock mattress settings, press and hold the LOCK button until the LED light turns off.

Note: MAX INFLATE mode is still functional in "LOCK" mode.



Quartet™ TA Quick Reference Guide



MAX INFLATE

Press the MAX INFLATE button and the mattress will inflate to its maximum volume within 60 seconds (MAX INFLATE automatically deactivates after 15 minutes).

Note: A series of beeps will sound every three minutes as a reminder that MAX INFLATE mode has been activated.



FIRMER



SOFTER

PATIENT COMFORT CONTROL LEVEL

The system is designed for patients weighing between 35-500 lbs., or up to 1000 lbs. on the bariatric model. Press the COMFORT CONTROL arrows to soften or firm the mattress to the patient's preferred comfort level.



SELECT

The SELECT button cycles between the *Pulse* and *Fowler* therapy modes. To select the desired mode, press the SELECT button until that mode is shown on the display. Use the ADJUST SETTINGS arrows to select pulsation frequency and duration options.

Note: For custom settings only, FOWLER mode will enable the mattress to transition to an upright position when the head of bed is articulated to 35° or greater.



ALARM SILENCE

Press the ALARM SILENCE button the control unit "beeps" to mute audio alarm.

BED EXIT ALARM

To enable the bed exit alarm, press and hold the ALARM SILENCE button until the control unit "beeps" and follow the instructions on the screen.



POWER STATUS

Notification that the unit is running off electrical power and the battery is charging.



BATTERY POWER

Notification that the unit is running on battery and the batteries are discharging. Displays percentage of battery capacity when the unit is running on either AC power or battery power.



LOW PRESSURE

The micro-controller will activate an audio/visual signal to alert the caregiver by flashing "LOW PRESSURE" on the digital display. Once the low-pressure problem is fixed, the audio/visual signal will cease.

SAFETY AIR PAD

The bottom safety air pad can be manually deflated by disconnecting the bottom pad deflate connector which is located at the bottom right corner (patient's right) of the mattress.

CPR FUNCTION

To deflate the mattress for a CPR procedure, disconnect the hose from the control unit.



Please read user manual for operating instructions and www.oncaremedical.com for additional product information.

For further questions, please call your local UHS office or call **1.800.847.7368**.