



EZ Way, Inc.

"Your Total Patient Handling Solution"

Smart Lift®

Mobile Full Body Patient Floor Lift in weight capacities of 500, 600, & 1,000 lbs.



EZ Way's proprietary circuit board system includes features to provide customers with added intelligence and safety features including:

Backlit LCD screen

During operation, EZ Way's LCD screen is illuminated for ease in viewing lift readout information.

Usage time captured

Lift usage time is captured in minutes, great for tracking frequency of use.

Overload protection

Should the actuator be overloaded due to lifting excess weight, lifting against an obstruction, or mechanical failure, "Actuator Overload" will appear on the LCD screen, the stand will cease lifting, and will only lower. A great feature to protect both patients and the lift.

Weight Overload

This feature will display "Weight Overload" on the LCD screen and will only allow the lift to lower.

Weight recall

On units with scale, the last weight taken is retained until a new weight is taken, allowing caregivers easy reference in busy working environments.


Smart Lift® Applications:

- Transfer patients to/from bed, commode, wheelchair, chair or from the floor.
- Reposition or turn patients in bed using specialty accessories.
- Weigh patients using optional scale.
- Elevate limbs using Limb Lifting Straps.



Models:

- L500PN** - 500 lb. Smart Lift®
- L500PS** - 500 lb. Smart Lift® w/scale
- L600PN** - 600 lb. Smart Lift®
- L600PS** - 600 lb. Smart Lift® w/scale
- L1000PN** - 1000 lb. Smart Lift®
- L1000PS** - 1000 lb. Smart Lift® w/scale
- LLB500PN** - 500 lb. Smart Lift® Low Profile Leg
- LLB500PS** - 500 lb. Smart Lift® Low Profile Leg w/scale
- LLB600PN** - 600 lb. Smart Lift® Low Profile Leg
- LLB600PS** - 600 lb. Smart Lift® Low Profile Leg w/scale

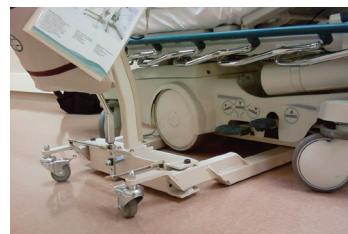

Low Profile Base

Optional low base for use with low clearance beds and stretchers.



Floor Lifting

Ability to lift patients from the floor.

"Your Total Patient Handling Solution"



DISPOSABLE SLING

Single patient use for infection control.

Available from Small to XXL, with and without head support.



WIPEABLE SLING

For infection control. Less frequent need to launder.

Available from Small to XXXL.



MULTI-PURPOSE SLING

For amputees or patients with larger lower extremities.

Available from Small to XXXL.

Additional Sling styles available.



CONTROL CONSOLE WITH LCD READOUT

Battery gauge, function indicator, dual up/down control, digital scale readout.

Smart Lifts have a
10 YEAR LIFE CYCLE

Average Annual Parts Cost:

\$38.22

Additional Features:

- 10 Year Frame / 3 Year Component Warranty
- Highest Weight Capacities in the Industry
- 3 Position Base Opening for Function Around Wheelchair, Chair and Commode
- Battery Powered; Always Available for Use
- Pneumatic Hand Control
- Digital Scale Available
- Low-base, 2", 4" and 5" Casters Available for Low Clearance of Special Beds or Ease of Function on Carpeted Flooring
- Redundant Controls
- Emergency Stop Safety Control
- Emergency Lowering Option
- Anti-pinch Feature
- 24-hour Service Hotline
- Includes 2 Batteries, Remote Charger and One Medium Deluxe Sling

	500 lb. Unit	600 lb. Unit	1000 lb. Unit
Battery Power	2 - 24 Volt lead acid		
Remote Charging Unit	1 - 120V UL listed		
Base Dimensions w/ 3" Standard Caster	24" w x 48" l x 4.75" h	26" w x 48" l x 4.75" h	33" w x 66" l x 4.75" h
Base Height w/ 2" Caster	3.25"		
Base Height w/ Low Base	2.25"		N/A
Sling With Unit	Size Medium Deluxe		Size XL Deluxe
Sling Bag	Standard		
Digital Scale	Optional		



EZ Way is proud to be a US manufacturer of Safe Patient Handling products



PO Box 89, 807 E. Main • Clarinda, IA 51632 • 1-800-627-8940

Form #2-238
Rev. 01/09/2020

WARNING: Cancer and Reproductive Harm - www.P65Warnings.ca.gov.