

# POWER VAC+

INTERMITTENT ASPIRATOR



## Power & intermittent suction in a small package!

### ***It's Powerful ...***

We've put power and higher flow into a small package. Our high performance 1/5 horsepower pump delivers a suction range of 0-200 mmHg and a full vacuum of 279 mmHg.

### ***It's Lightweight & Compact...***

An average of 42% lighter & 44.8% smaller than comparable units, it fits into areas where space is at a minimum and is easy to carry.

### ***It's Convenient ...***

The unique canister holder is moveable, allowing the caregiver to rotate the canister to either side of the unit for best placement near the point of use.

### ***It's Durable ...***

The ABS plastic cabinet is durable, attractive, easy-to-clean and infection control friendly.

### ***It's Protected ...***

A dual protection system consists of a float shut-off in the canister lid and an Inline Hydrophobic Filter to prevent back flow of fluid as well as airborne contaminants from entering the pump.


### ***It's Quieter ...***

Our unique design allows a quieter operation for around the clock use.

### ***It's Field Serviceable & Dependable ...***

The PowerVac+ Intermittent Aspirator has a unique patented modular designed regulator. No need to replace the full unit, only the problem module.



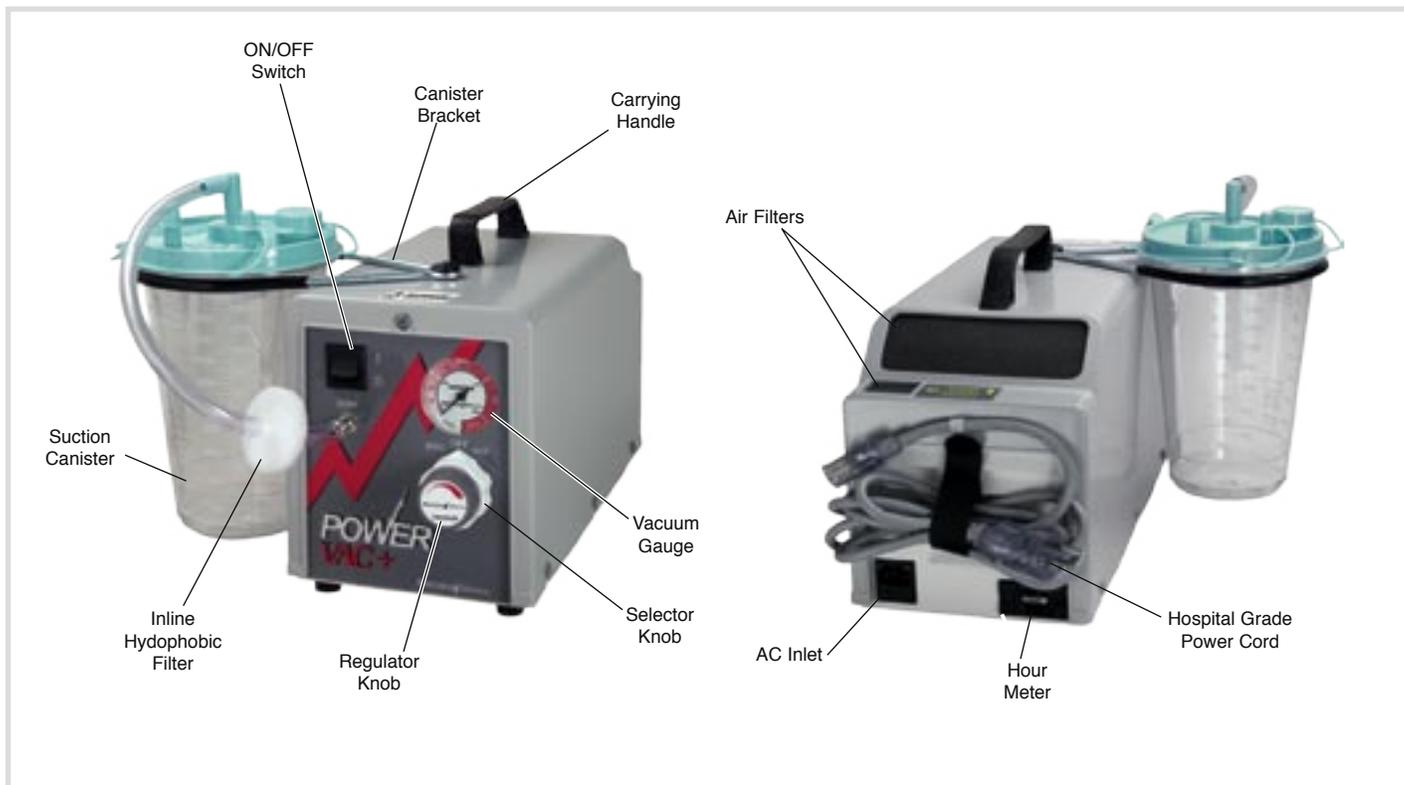

## SPECIFICATIONS

<b>Dimensions:</b>	Length 15 <sup>7</sup> / <sub>8</sub> " (40.3 cm), Width 7 <sup>3</sup> / <sub>8</sub> " (18.7 cm), Height 10 <sup>1</sup> / <sub>8</sub> " (25.7 cm)
<b>Weight:</b>	21.7 lbs. (9.8 kg)
<b>Shipping Weight:</b>	26.1 lbs. (11.8 kg)
<b>Motor:</b>	120 VAC, 60Hz, 5 Amps

<b>Vacuum Range:</b>	0-200 mmHg with Full Vacuum of 279 mmHg
<b>Horsepower:</b>	1/5
<b>Ambient Temperature Range:</b>	41°F to 95°F (5°C to 35°C)
<b>Storage Conditions:</b>	Temperature: -4°F to 104°F (-20°C to 40°C); Humidity: Max 95% non-condensing

## ORDERING INFORMATION

Part#	Description
PM63	Power Vac+ Aspirator with Hospital Grade Power Cord & 2000cc canister
1955	Patient Suction Tubing (6ft.)
502540	Canister Bracket

Part#	Description
502687-12	Canister (with Lid), 2000cc (Case of 12)
502690	In-line Hydrophobic Filter & Tube Assembly
PM60-1	Roll Stand
PM60-2	Roll Stand with Suction Canister Bracket



ISO 9001  
Certified

Precision Medical, Inc.  
300 Held Drive • Northampton, PA 18067  
Phone: 610-262-6090 Toll-Free Phone: 1-800-272-7285  
Fax: 610-262-6080 Toll-Free Fax: 1-800-352-1240  
Visit us on the Web: [www.precisionmedical.com](http://www.precisionmedical.com)