


ELITE ACUTE CARE BED WITH POWER ASSIST DRIVE

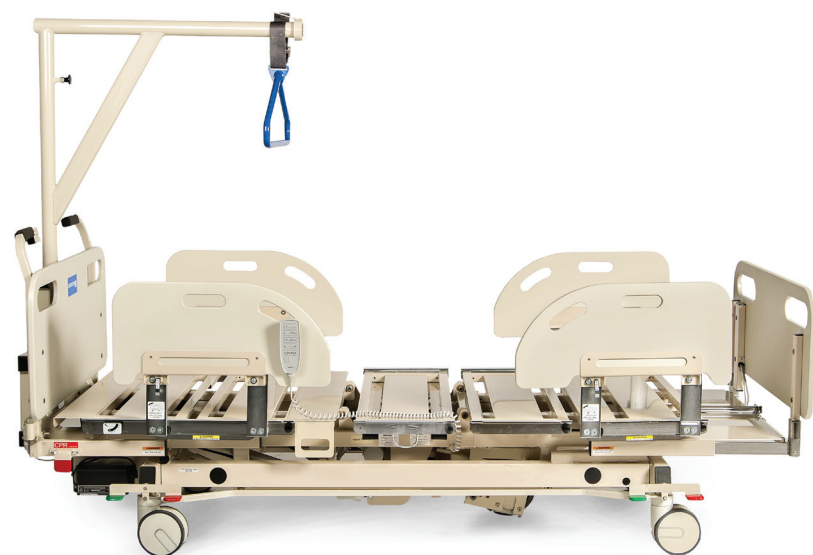

- PATIENT WEIGHT CAPACITY:** 1,000 pounds
- FOWLER BACK REST ANGLE:** 0-65 degrees
- FOLEY BAG HOLDER:** One attachment each side of bed
- RESTRAINT ATTACHMENT:** One attaching point, each side of bed
- POWER ASSIST DRIVE:** Engage P.A.D. by lowering bed deck to lowest position. Disengage by raising bed deck
- VARIABLE SPEED CONTROL:** Switch embedded in right side head end. Push handle for forward/reverse Power Assist Drive functions
- TRENDELENBURG/REVERSE TRENDELENBURG:** 9 degrees
- KNEE FLEX/GATCH:** 0-20 degrees
- DECK LOW POSITION:** 17"
- DECK HIGH POSITION:** 26"
- ADJUSTABLE DECK LENGTH:** 84"-92"
- ADJUSTABLE DECK WIDTH:** 39"-48"-54"
- SIDE RAILS (4):** 17" high above deck
- CPR FOWLER BACKREST RELEASE:** Head end, left and right
- FOUR WHEEL LOCK AND STEER:** Engage wheel locks or steer at any caster wheel
- INTEGRATED PUSH HANDLES:** Two, at head end
- PATIENT HELPER/TRAPEZE:** Swings 90 degrees from center
- CONTROL PENDANT:** Illuminated, five-function, self-storing
- CAREGIVER CONTROL:** Controls all functions and locks out pendant functions
- BED POSITIONING/ARTICULATION:** High/low, fowler backrest, Trendelenburg/Reverse Trendelenburg, knee flex/leg lift, cardio chair
- ANGLE INDICATOR:** Two at head end, one left, one right
- UNDER BED ILLUMINATION:** LED illumination
- OUT OF BED EXIT ALARM:** Sounds at bed with option to sound at nurse call station using DB 37 cable or 1/4" phone connector
- INTEGRAL SCALE WEIGHING SYSTEM:**
- Displays accurate weight in any articulated position including Trendelenburg and Reverse Trendelenburg and in any caster orientation
 - One convenient control box contains both the scale controls and the caregiver key pad
 - Scale display rotates - HIPAA compliant
 - Powered by bed electrical system and battery backup. Does not require batteries and will function when disconnected from A/C power
 - Weight change mode
 - Dual weight display
 - Freeze mode
 - Backlit display
 - Accurate to within 0.2% or 2 pounds, whichever is greater

PATENT PENDING
DESIGNED and MANUFACTURED IN USA

- WHY BUY GENDRON:**
- 25 years experience serving the bariatric patient
 - Designed and manufactured in the USA
 - Complete bariatric product portfolio

GENDRON


800.537.2521 • sales@gendroninc.com
www.gendroninc.com

GENDRON

TOTAL
SOLUTIONS
FOR
BARIATRIC
PATIENT
CARE


MOBILITY
SAFETY
COMFORT
DIGNITY

The Elite Acute Care Bed with Power Assist provides your bariatric patients with all of the above in an easy-to-use design, engineered to allow you to give your patients the best possible care. Now, with the addition of Power Assist Drive, the Elite bed provides maximum mobility and ease of use for staff and transport team. Easily navigate and maneuver with the Power Assist Drive.

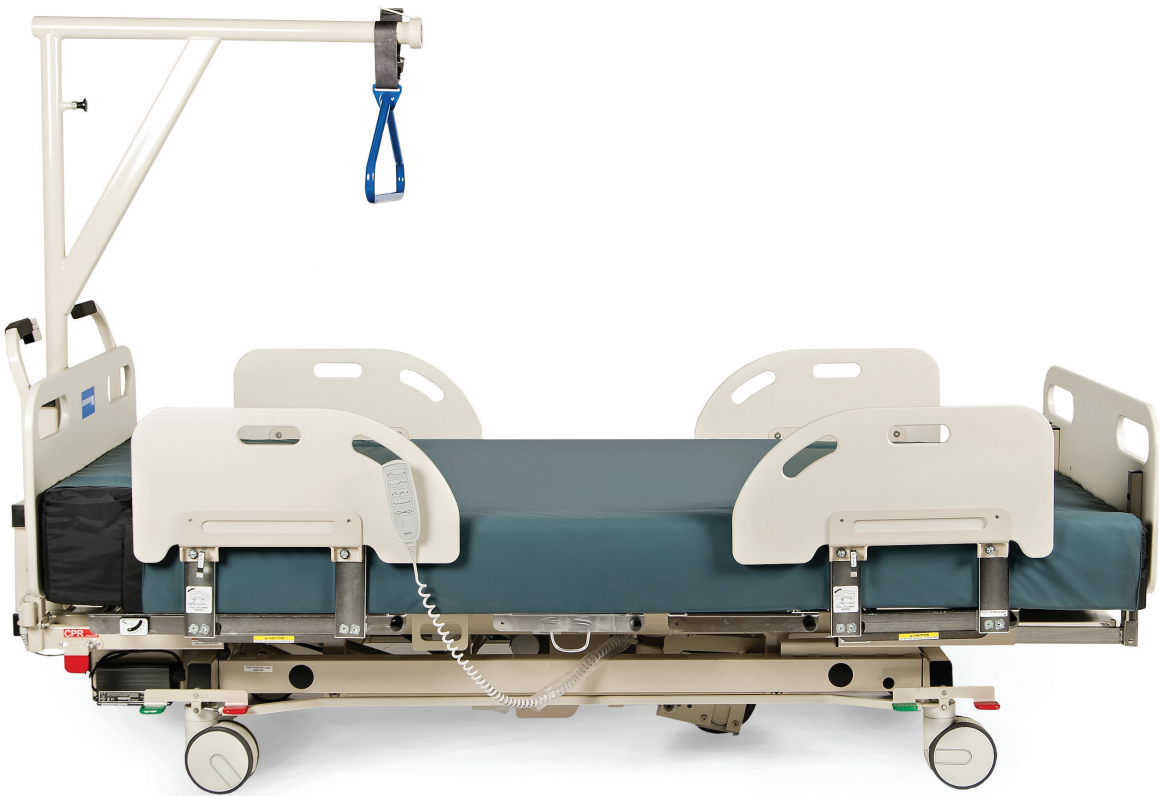

ELITE BARIATRIC ACUTE CARE BED
FEATURING POWER ASSIST DRIVE

Elite Acute Care Bed with Power Assist Drive Features




VARIABLE SPEED SWITCH, FORWARD AND REVERSE

ON/OFF POWER SWITCH FOR POWER ASSIST DRIVE

NURSE CALL CONNECTIONS

QUICK ELECTRICAL DISCONNECT FOR HEADBOARD REMOVAL

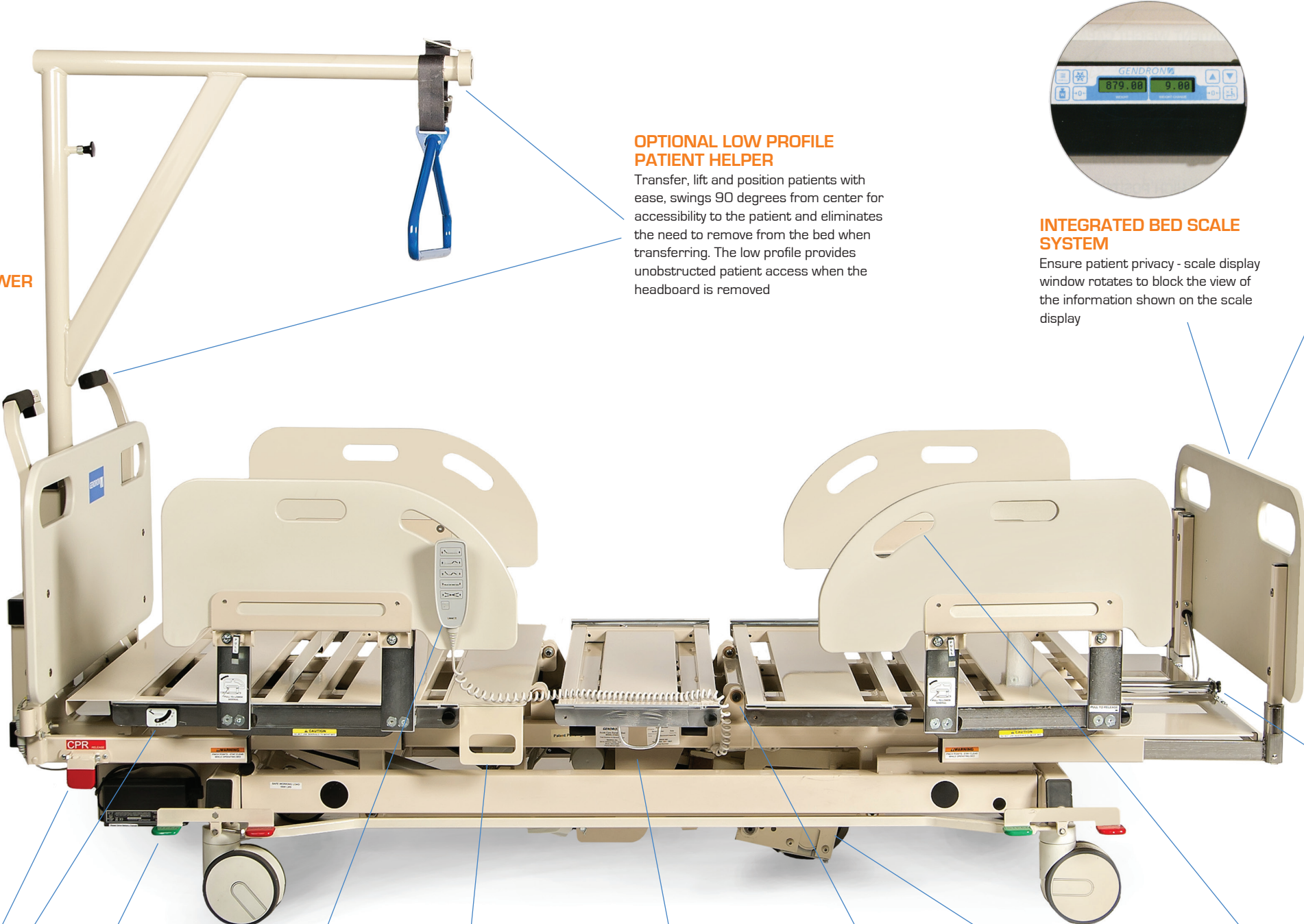

OPTIONAL LOW PROFILE PATIENT HELPER
Transfer, lift and position patients with ease, swings 90 degrees from center for accessibility to the patient and eliminates the need to remove from the bed when transferring. The low profile provides unobstructed patient access when the headboard is removed


INTEGRATED BED SCALE SYSTEM
Ensure patient privacy - scale display window rotates to block the view of the information shown on the scale display

CAREGIVER CONTROL
Located at foot end of bed. Any or all hand pendant functions can be locked-out at the caregiver control to restrict patient access


FEATURES INCLUDE:

- BED POSITIONING/ ARTICULATION**
Provide proper clinical positioning from Trendelenburg and Reverse Trendelenburg to knee flex and cardio chair. Easy-to-use button controls; one button cardio chair positioning and one button return to level positioning
- OUT OF BED EXIT ALARM**
Alerts you when patient attempts to exit bed. Programmable alarm sounds at the bed, with an option to connect to the nurse station


ADJUSTABLE DECK LENGTH
Ensure the comfort of your tall patients - one touch extension from 84" to 92" long

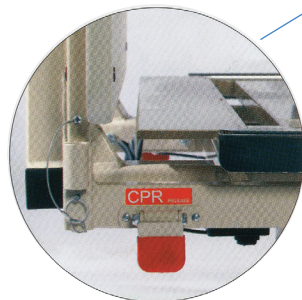


CPR FOWLER BACKREST RELEASE
Respond faster to patient needs: lower backrest to flat within six seconds for CPR




FOWLER BACKREST
Minimum 65 degree angle in raised position

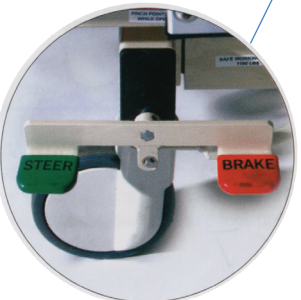


FOUR WHEEL LOCK AND STEER
Steer safely and precisely - engage wheel lock and steer at any caster wheel




CONTROL PENDANT
Adjust settings confidently even in dark rooms. Illuminated controls; five-function, self-storing at head end, left and right side

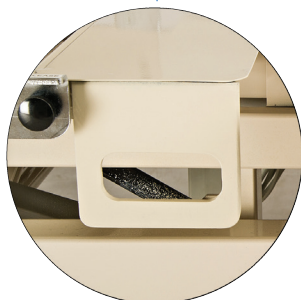


WRIST RESTRAINT ATTACHMENT POINT
Located on both sides of the bed




FOLEY BAG HOLDER
Located on both sides of the bed

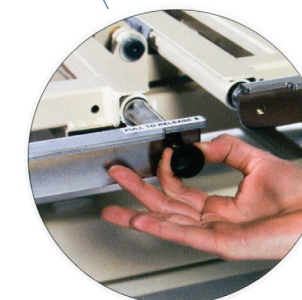

ADJUSTABLE DECK WIDTH
Transport safely and provide comfort - one touch extension from 39" to 48" to 54" wide



POWER ASSIST DRIVE MECHANISM
Incorporates mechanical electric brake manual release



DUAL HAND HOLDS
On all four side rails for patient safety