

Cub® Pediatric Crib

Medical

Rental Solutions Delivered by



1 800 STRYKER Option 3

Universal | Hospital | Services®



Easily-lowered siderails for unobstructed access

Swing-open access doors allow quiet, convenient patient care

Retracting backrest helps eliminate head-end gap


Integrated scale system obtains accurate weights with the child in any position

Retractable fifth wheel helps improve mobility



Optional scale system provides accurate weights with the touch of a button.

Accent Color Options

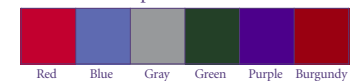

Safety for Smaller Patients

Responsible pediatric care demands a stretcher specifically designed around safety, access and ease of use. The Cub pediatric crib supports safe and efficient care practices while the warm, friendly appearance cultivates a healing environment for small children and their families.



Optional protective top equipped with intermediate support position on each side.



Retracting backrest reduces the gap against the siderail while swing-open access doors allow easy patient care.



Lowered siderails provide unobstructed access.



190 Hydraulic Model



190 Fixed Height Model

Medical

3800 E. Centre Ave.
Portage, MI 49002
t: 269 329 2100 f: 269 329 2311
toll free: 800 787 9537

www.stryker.com

Standard Features

- 3" (8 cm) vinyl mattress
- Manually-operated retracting backrest
- Manual foot prop
- Siderail with 9" (23 cm) mandatory stop
- Counterbalance system on siderail for easy raising and lowering
- Oxygen bottle holder
- Fixed foot-end rail

Hydraulic Model: All Standard Features Plus

- 6" (15 cm) Omni Surface™ casters
- Central brake control
- Dual side-mounted foot control hydraulics with uni-lower pedal
- Hands-free Trend./reverse Trend.

Specifications

	190 Hydraulic		190 Fixed Height	
• Model Number				
• Overall Length & Width				
No Accessory Supports	65.5" x 37"	(166 x 94 cm)	65.5" x 37"	(166 x 94 cm)
With Accessory Supports	72.25" x 37.5"	(184 x 95 cm)	72.25" x 37.5"	(184 x 95 cm)
• Weight Capacity	400 lbs	(181 kg)	400 lbs	(181 kg)
• Height Range (to litter top)				
High	40"	(102 cm)	33"	(84 cm)
Low	33"	(84 cm)	33"	(84 cm)
• Litter Positioning				
Manual Activated	0° / 25° / 40° / 55°		0° / 25° / 40° / 55°	
Pneumatic Assist	0° - 50°		0° - 50°	
Trend./Reverse Trend.	±9.5°		0° / 6° / 12° with manual foot prop	
Foot Section Angle	0° / 6° / 12°		0° / 6° / 12°	
• Patient Surface	29.5" x 57"	(75 x 145 cm)	29.5" x 57"	(75 x 145 cm)
• Caster Diameter	6"	(15 cm)	5"	(13 cm)

Stryker reserves the right to change specifications without notice.

Warranty

- One year parts, labor and travel or two years parts

Optional Features

- Integrated scale system
- 4" (10 cm) Ultra Comfort Mattress
- Siderail without 9" (23 cm) mandatory stop
- Access doors at all four corners
- Fully-lowerable foot-end rail – provides unobstructed access around entire unit when all siderails are down
- Accessory supports with roller bumpers (Head, foot or both ends)
- Chart shelf
- Retractable fifth-wheel steering
- IV caddy
- Monitor tray (Not compatible at foot end when scale option is selected)
- Oxygen bottle retaining ring
- IV poles
 - Permanent folding two- or three-stage
 - Standard removable
- Pneumatic assist backrest
- Retractable protective top
- Retractable protective top with intermediate support position
- Upright oxygen bottle holder
- Siderail pads

